- 1. The sheild consists of two pieces:
- The primary shield




- And the headband




2. Start by bending the top part of the shield and crease to form a seam



3. Insert the tabs on the shield into the slots.

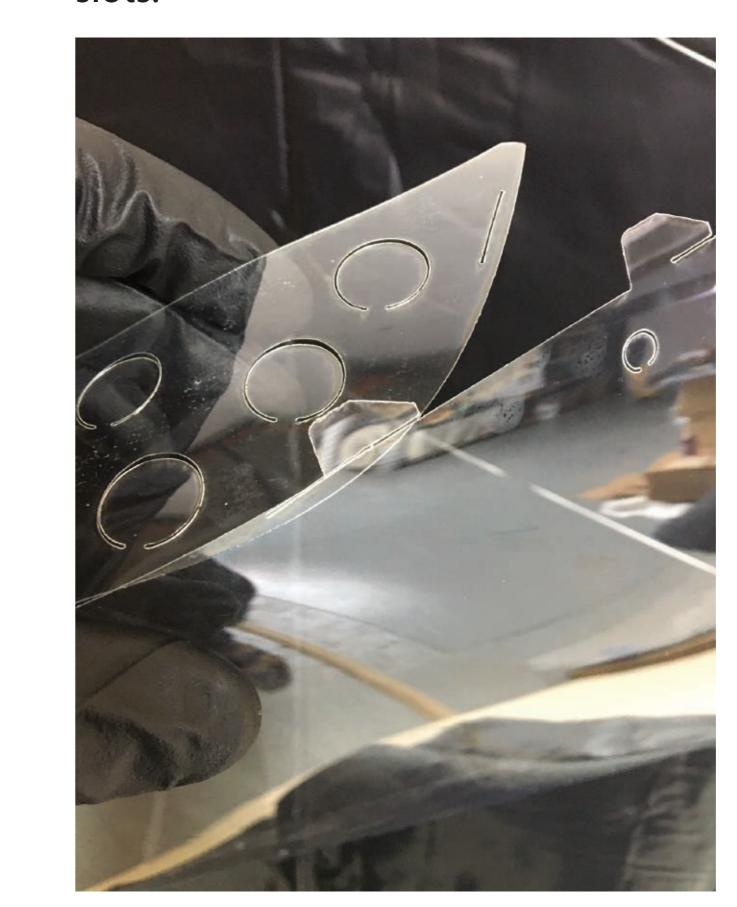


4. Bend the tabs forward and crease firmly to lock the tabs in place.

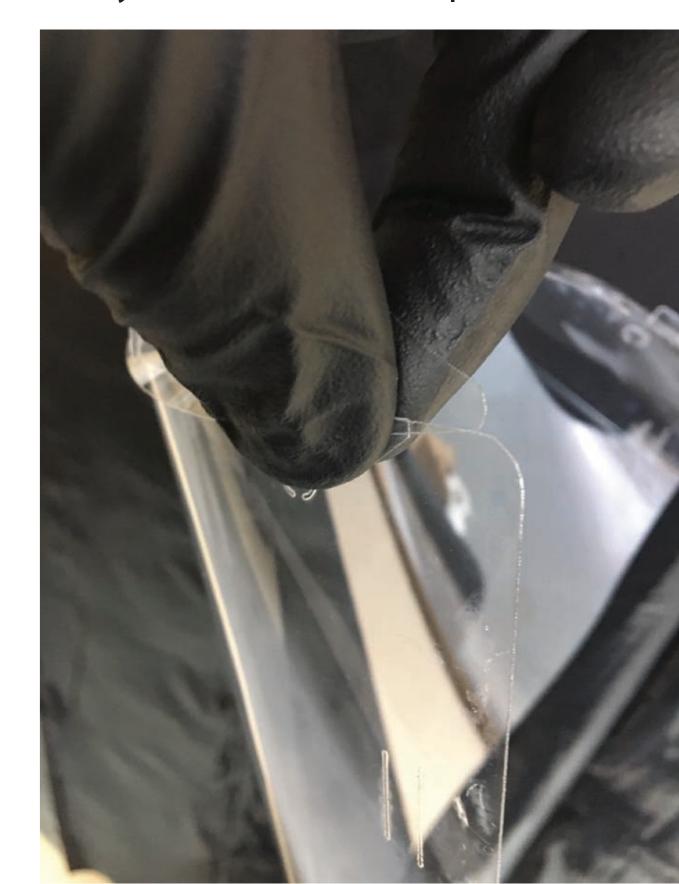

5. The main shield is now assembled




5. Push open and crease the vent holes to the desired position.

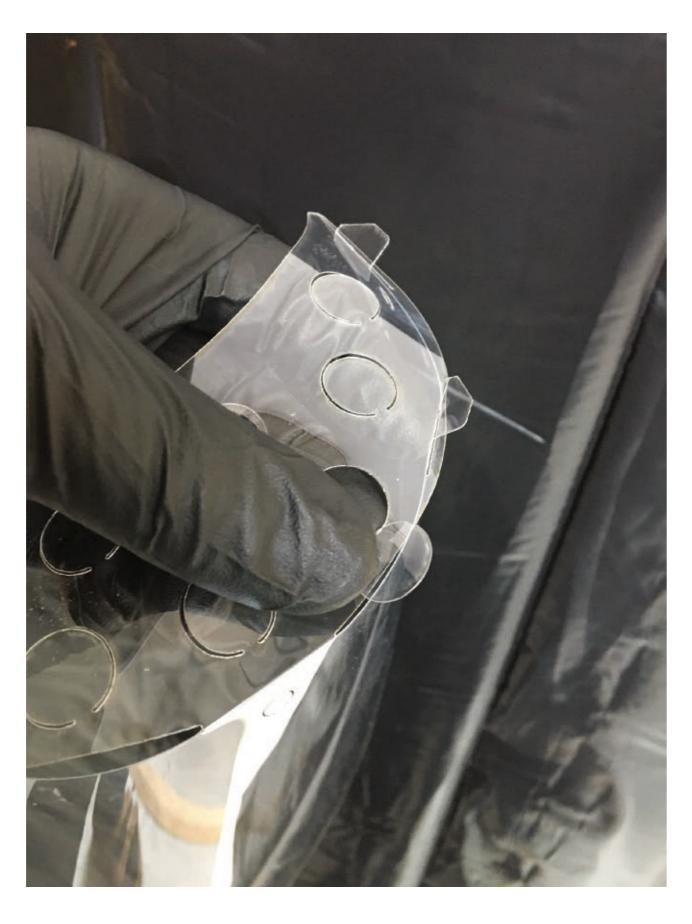


7. Push a tab of the headband through the outer slot on the side of the shield



8. Bend the tab back toward the front of the shield and crease sharply

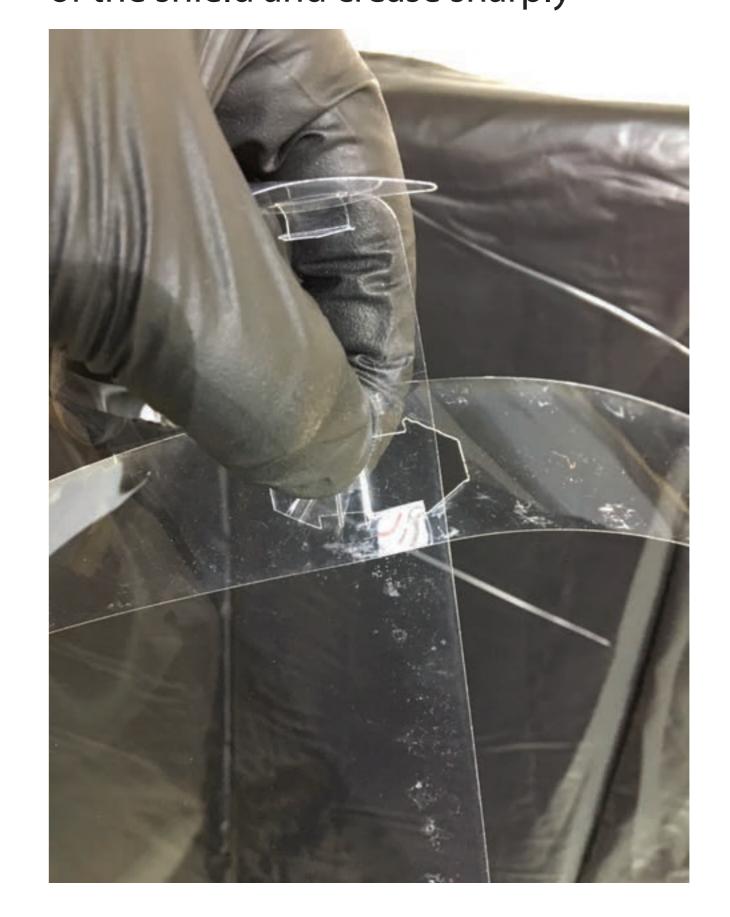

9. Insert the tab through the inner slot on the shield. (Alternatively, this connection can be stapled in place)

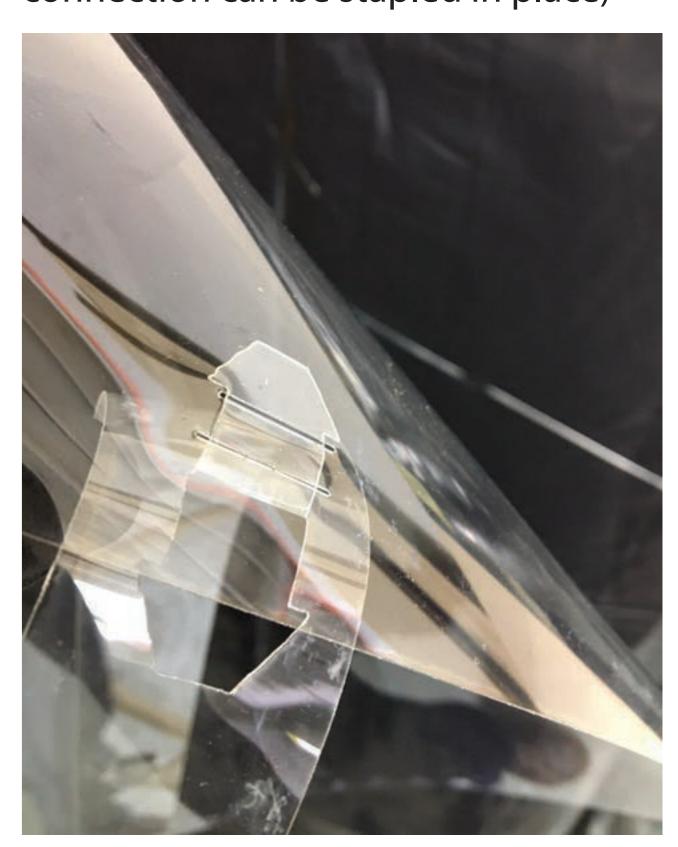

PLEASE SANITIZE APPROPRIATELY BEFORE USE